



Surgical Accessories

Unomedical Pty Ltd
11-17 Wilmette Place,
Mona Vale, NSW, Australia 2103
T (02) 9997 8033
F (02) 9999 4031

New Zealand
Customer Service: 0800 161 777
Product Support: 0800 445 400
www.unomedical.com

Local response, global expertise

Unomedical leadership


Utility, comfort and performance are keynote characteristics of the high quality yet value-for-money Unomedical range of Surgical Accessories.

Drawn from the very best available both here and internationally, the Unomedical range is developed specifically for Australasian hospitals.



Surgical Accessories Products

- Drapes
- Anti-Fog Wipes
- Tourniquet Tubes
- Camera Covers
- Laparoscopic Insufflation Tubing
- Maxxim Brush

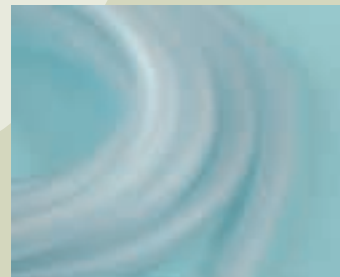

Drapes

Drapes

Polyethylene, clear. Folded so as to allow one handed opening.
Peel pack, sterile.

REF.	Drape Size	Qty / Box
298.48.061	48 cm x 61 cm (18" x 24")	200
298.61.092	61 cm x 92 cm (24" x 36")	200
298.71.102	71 cm x 102 cm (28" x 40")	200
298.92.120	90 cm x 120 cm (35" x 47")	200
298.99.999	122 cm x 137 cm (48" x 54")	100

Mayo Table Sleeve

Polyethylene,sealed at one end, clear: Folded, so as to allow one handed opening.
Peel pack, sterile.

REF.	Drape Size	Qty / Box
300.61.092	61 cm x 92 cm (24" x 36")	100

Polysheet

Polyethylene, blue. Folded so as to allow one handed opening.
Supplied bulk and clean.

REF.	Drape Size	Qty / Box
363.91.150	91 cm x 150 cm (36" x 59")	250



Other Surgical Accessories



Anti-Fog Wipes

REF.	Description	Qty / Box
RT-100	Reso Clear	24



Tourniquet Tubes

Tourniquet tubes can be used as leverage in tying off major blood vessels. Used in cardiothoracic surgery. Sterile, peel pack.

TOURNIQUET TUBE.

REF.	OD	Length	Qty / Box
605.16.220	5.4mm	22cm	100

ASD TOURNIQUET TUBE.

REF.	OD	Length	Qty / Box
605.24.180	8.0mm	18cm	100

Camera Covers



REF.	Description	Qty / Box
2904	Camera Cover Universal	25
3501	Camera Cover Closed	25



Laparoscopic Insufflation Tubing with Filter

Prevents cross contamination between patient and equipment. Latex free

Filter media is a barrier against

- Tuberculosis
- Hepatitis A
- Hepatitis B
- Hepatitis C

Unomedical Insufflation Tubing System, manufactured in Australia with significant user advantages:

Hydrophobic filter – for performance and safety

- Low deadspace bacterial and viral filter offering optimum hydrophobic filtration performance
- HEPA rated hydrophobic media allows effective removal of airborne and waterborne particulate, viral and bacteria
- Overall design maximises air/gas flow rates and lowers pressure drop
- Bi-directional filter media eliminates orientation concerns
- Lightweight housing ensures tubing integrity
- Filter media has a filter efficiency of 99.999%

Special Soft Tubing – for performance

- Tubing manufactured from high quality plastic
- Tubing properties similar to silicone tubing, soft and flexible, which prevents coiling and kinking

Technical Data

Filter Media	Filter media meets ULPA standards for air/gas delivery with 99.999% filter efficiency at test flows of 40 L/min
Pressure Hold	1 psi (0.069 bar) hold for 10 minutes
Minimum Air Flow Rate (Lpm at 1 psi)	≥ 15
Effective Filtration Area	18.5cm ²
Sterilisation	ETO

REF.	Description	Tubing Length	Qty / Box
515.10.300	Insufflation tubing with Silicone Connector	3 metres	20
515.20.300	Insufflation tubing with filter 2 male Luer Locks	3 metres	20
515.30.300	Insufflation tubing with filter and Universal Connector	3 metres	20
515.50.300	Insufflation tubing with filter and 15/22mm Connector	3 metres	20



Maxxim Brush

INFECTION PREVENTION WITH A SOFT TOUCH

Good hygiene forms the basis for preventing infection in the hospital. with good reason, since the majority of hospital infections are the result of insufficient hand and skin hygiene. The introduction of the Maxxim Brush in your hospital will enable you to establish the right conditions for optimal preoperative hygiene of hand and skin of staff and patients.

A separate nail cleaner is also enclosed, the fact that brush, soap and nail cleaner are combined in one package means that the Maxxim Brush also reduces handling costs. All these features combine to make Maxxim Brush the complete solution for good preoperative hygiene of hand and skin of staff and patients.

The Maxxim Brush comes in a well sealed package which keeps the soap stable but remains very simple and hygienic to open.

The handy dispenser box contains 25 brushes and can be attached to the wall beside the wash basin and/or in the OR.

A light weight aluminium holder is also available. The strictest of quality controls (ISO 9001) are applied throughout the production process and the Maxxim Brush complies fully with the very latest standards.



The disposable Maxxim Brush has many advantages which traditional reusable brushes simply cannot offer. The Maxxim Brush is very convenient to use. Its soft bristles for nails, knuckles and skin clean effectively while being kind to the skin. The attached sponge gently cleans the hands, forearms and skin avoiding the risks and discomfort of skin irritation and abrasion.

Physical and chemical details:

	Chlorhexidine	PVP iodine
pH value	± 6,0	± 5,35
density (g/cm ³)	± 1,02	± 1,07
appearance	amber liquid	brown liquid



All the benefits of the disposable Maxxim Brush at a glance:

- significant reduction of microorganisms
- deeper-lying resident bacterial flora are not removed unnecessarily
- causes no microtraumas on the skin, reducing risk of infection
- soft and kind to the skin, encouraging more thorough use
- effective, sensitive and conveniently shaped brush, comprising two layers of ultrasoft material
- balanced supply of soap through the squeeze-through baseplate
- no waste of soap
- no need for resterilization procedures, saving time, work and money.

There are three types of Maxxim Brush:

- 1 DRY STERILE**
 - sterilized using ethylene-oxide (ETO)
 - without soap
 - easy to open
- 2 IMPREGNATED WITH CHLORHEXIDINE GLUCONATE IN DETERGENT**
 - contains 4% chlorhexidine digluconate
 - broad spectrum of action
 - good residual working
 - mild cleansing agent: minimal chance of skin irritations
 - enough soap for 3-5 minutes' washing
- 3 IMPREGNATED WITH PVP IODINE IN DETERGENT**
 - contains a minimum of 1% iodophor
 - broad spectrum of action
 - fast and efficient action
 - mild cleansing agent: minimal chance of skin irritations
 - enough soap for 3-5 minutes' washing

REF.	Brush Type	Qty per dispenser box	Boxes per case
255A480	Dry sterile	25	15
255A481	Impregnated with chlorhexidine gluconate in detergent	25	15
255A482	Impregnated with Povidine iodine in detergent	25	15
04823	Scrub brush holder	-	1